
# FORA<sup>®</sup> Premium V10

Blood Glucose Monitoring System

*Sistema de vigilancia de glucose sanguinea*

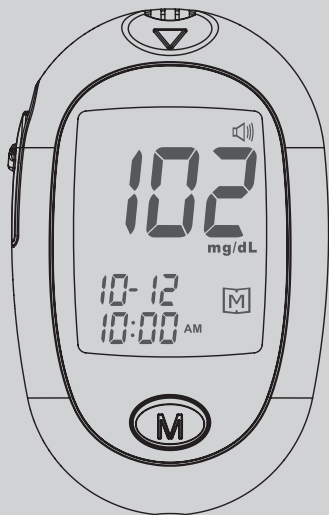**Owner's Manual**  
***Manual del Propietario***

# Dear FORA Premium V10 System Owner:

Thank you for purchasing the **FORA Premium V10** Blood Glucose Monitoring System. This manual provides important information to help you use the system properly. Before using this product, please read and follow the instructions in this Owner's Manual and any inserts for the test strips and control solution.

Regular monitoring of your blood glucose levels can help you and your doctor gain better control of your diabetes. Due to its compact size and easy operation, you can use the **FORA Premium V10** Blood Glucose Monitoring System to easily monitor your blood glucose levels.

If you have questions, just call ForaCare's Customer Care Service Center toll-free at 1-888-307-8188 (8:30am-5:00pm PST, Mon.-Fri.). You can also visit [www.foracare.com](http://www.foracare.com).

# IMPORTANT SAFETY PRECAUTIONS

## READ BEFORE USE

- The meter and lancing device are for single patient use. Do not share them with anyone including other family members or use on multiple patients!
- All parts of the kit are considered biohazardous and can potentially transmit infectious diseases, even after you have cleaned and disinfected the components in the kit.

For more information, please visit

1. "FDA Public Health Notification: Use of Fingerstick Devices on More than One Person Poses Risk for Transmitting Bloodborne Pathogens: Initial Communication" (2010)

<http://wayback.archive-it.org/7993/201701111013014/http://www.fda.gov/MedicalDevices/Safety/AlertsandNotices/ucm224025.htm>

2. "CDC Clinical Reminder: Use of Fingerstick Devices on More than One Person Poses Risk for Transmitting Bloodborne Pathogens" (2010)

<http://www.cdc.gov/injectionsafety/Fingerstick-DevicesBGM.html>

1. Use this device **ONLY** for the intended use described.
2. Do **NOT** use accessories which are not specified by the manufacturer.
3. Do **NOT** use the device if it is not working properly or if it is damaged.
4. Do **NOT** use the equipment where aerosol sprays are being used or where oxygen is being administered.
5. This device does **NOT** serve as a cure for any symptoms or diseases. The data measured is for reference only.
6. Before using this device to test your blood glucose, read all of the instructions thoroughly, practice the test with the meter and perform all quality control checks as directed.
7. Please use this device in a dry environment, especially when synthetic materials are present ( clothing, carpets etc.). These materials give off static discharges, which may cause erroneous results.
8. Do **NOT** use this instrument in close proximity to sources of strong electromagnetic radiation, as these may interfere with the accurate operation.

9. Do not under any circumstances use the device on neonates or infants.
10. Keep the device and testing equipment away from young children. Small items such as the battery cover, batteries, test strips, lancets and vial caps are choking hazards.
11. Proper maintenance and periodically checking the meters quality through control solution tests, are essential to the longevity of your device. If you are concerned about the accuracy of your measurement, please contact your local customer service or place of purchase for help.

**KEEP THESE INSTRUCTIONS IN A SAFE PLACE**

# TABLE OF CONTENTS

<b>BEFORE YOU BEGIN</b>	<b>07</b>
Important Information	07
Intended Use	08
Test Principle	08
Contents of System	09
Meter Overview	10
Display Screen	11
Test Strip	12
Setting the Meter	13
<b>THE FOUR MEASURING MODES</b>	<b>16</b>
Control Solution Testing	17
Performing a Control Solution Test	17
<b>TESTING WITH BLOOD SAMPLE</b>	<b>20</b>
Preparing the Puncture Site	20
Alternative Site Testing	22
Performing a Blood Glucose Test	23
<b>METER MEMORY</b>	<b>26</b>
Reviewing Test Results	26
Reviewing Blood Glucose Day Average Results	27
Downloading Results onto a Computer	28
<b>MAINTENANCE</b>	<b>29</b>
Battery	29
Caring for Your Meter	30
Caring for Your Test Strips	33
Important Control Solution Information	34
<b>SYSTEM TROUBLESHOOTING</b>	<b>35</b>
Result Readings	35
Error Messages	36
Troubleshooting	38
<b>DETAILED INFORMATION</b>	<b>39</b>
<b>SPECIFICATIONS</b>	<b>40</b>

# BEFORE YOU BEGIN

## Important Information

- Severe dehydration and excessive water loss may cause readings which are lower than actual values. If you believe you are suffering from severe dehydration, consult a healthcare professional immediately.
- If your blood glucose results are lower or higher than usual, and you do not have any symptoms of illness, repeat the test. If you have symptoms or continue to get results which are higher or lower than usual, consult with your healthcare professional.
- Use only fresh whole blood samples to test your blood glucose. Using other substances will lead to incorrect results.
- If you are experiencing symptoms that are inconsistent with your blood glucose test results and you have followed all the instructions given in this owner's manual, contact your healthcare professional.
- We do not recommend using this product on severely hypotensive individuals or patients in shock. Readings lower than actual values may occur for individuals experiencing a hyperglycemic-hyperosmolar state, with or without ketosis. Please consult a healthcare professional before use.
- Critically ill patients should not be tested with a blood glucose monitoring system.

## Intended Use

The **FORA Premium V10** Blood Glucose Monitoring System is intended for use in the quantitative measurement of glucose in fresh capillary whole blood from your finger and the following alternative sites: palm, forearm and upper arm. It is intended for use by people with diabetes mellitus at home as an aid in monitoring the effectiveness of diabetes control program. It is not intended for the diagnosis of or screening for diabetes mellitus, and is not intended for use on neonates.

It is intended to be used by a single person and should not be shared.

The alternative site testing in the **FORA Premium V10** Blood Glucose Monitoring System can be used only during steady-state blood glucose conditions.

This meter has some speaking functions but is not intended for use by the visually impaired.

## Test Principle

Your system measures the amount of sugar (glucose) in whole blood. The glucose testing is based on the measurement of electrical current generated by the reaction of glucose with the reagent of the strip. The meter measures the current, calculates the blood glucose level, and displays the result. The strength of the current produced by the reaction depends on the amount of glucose in the blood sample.


This system uses enzyme of glucose oxidase (GOD) from *A. niger*.

# Contents of System

The **FORA Premium V10** system kit includes:

- 1 **Meter**
- 2 **Owner's Manual**
- 3 **Quick Start User Guide / Warranty Card**
- 4 **Protective Wallet**
- 5 **Daily Log Book**
- 6 **Two AAA Batteries**

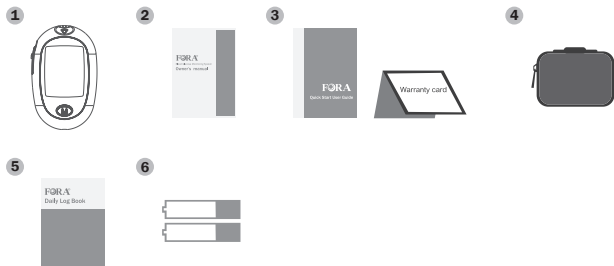
Test strips, control solution, sterile lancets and lancing device may be purchased separately. Please make sure you have those items needed for a blood glucose test beforehand.

If you wish to transfer data from your meter to a computer (see "Downloading Results onto your computer" p.28), you will need to purchase an RS232 cable, this cable is not provided in the kit. Please contact your local customer service or place of purchase to obtain this cable.

## NOTE:

If any items are missing from your kit or opened prior to use, or to obtain test strips and control solution, please contact the place of purchase or call Customer Service for assistance.

## Meter Overview

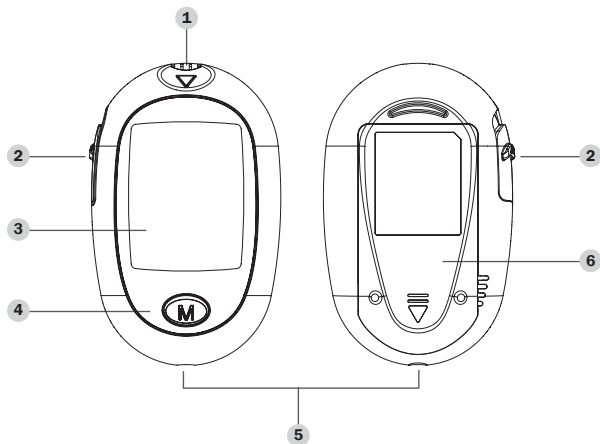**1 Test Strip Slot**

Insert test strip here to turn the meter on for testing.

**2 Test Strip Ejector**

Eject the used strip by sliding up this button.

**3 Display Screen**

**4 M Button**

Enter the meter memory and silence a reminder alarm.


**5 Data Port**

Download test results with a cable connection.

**6 SET Button**


Enter SET mode and confirm the meter settings.

## Display Screen

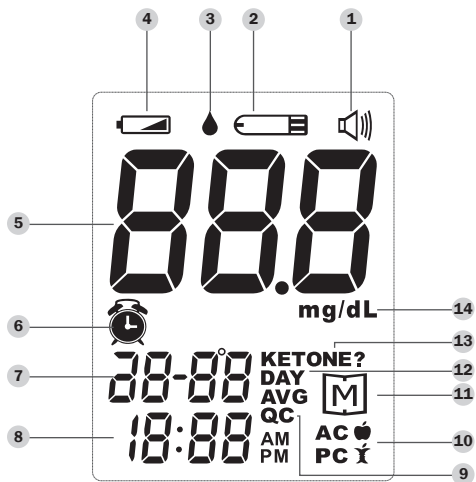- |   |                    |    |                       |
|---|--------------------|----|-----------------------|
| 1 | Volume Symbol      | 8  | Time                  |
| 2 | Test Strip Symbol  | 9  | Control Solution Mode |
| 3 | Blood Drop Symbol  | 10 | Measuring Mode        |
| 4 | Low Battery Symbol | 11 | Memory Symbol         |
| 5 | Test Result        | 12 | Day Average           |
| 6 | Alarm Symbol       | 13 | Ketone Warning        |
| 7 | Date               | 14 | Measurement Unit      |

## Test Strip
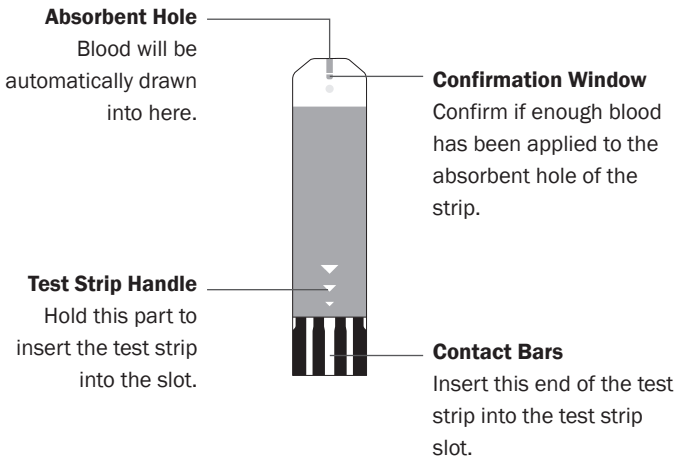
### ATTENTION:

**The front side of test strip should face up when inserting test strip.**


Test results might be wrong if the contact bar is not fully inserted into the test slot.

### NOTE:

The **FORA Premium V10** monitor should only be used with **FORA Premium V10** Test Strips. Using other test strips with this meter can produce inaccurate results.

## Setting the Meter

Before using your meter for the first time or if you change the meter battery, you should check and update the settings. Make sure you complete the steps below to change your settings.
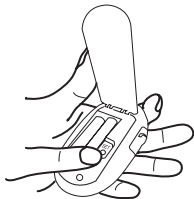
### Entering the Setting Mode


Start with the meter off (no test strip inserted). **Press SET.**

#### 1. Setting the date

With the year flashing, press **M** until the correct year appears. Press SET.

With the month flashing, press **M** until the correct month appears. Press SET.

With the day flashing, press **M** until the correct day appears. Press SET.
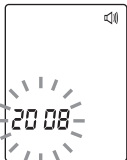
#### 2. Setting the time

With the hour flashing, press **M** until the correct hour appears. Press SET.


With the minute flashing, press **M** until the correct minute appears. Press SET.

### 3. Deleting the memory


With “dEL” and “M” on the display, press **M** and select “no” to keep the results in memory then press SET to skip.


To delete all the results, press **M** and select “yes” to delete all the memory records.
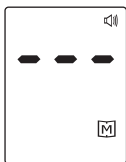
### 4. Setting the speaking volume

There are seven (7) speaking volume levels to choose from. Press **M** until the desired speaking volume appears. To confirm your selection, press SET.

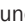Volume 0 indicates that the speaking function is turned off, and “” will not display during testing. Volume 1 to 7 indicates speaking volume from low to high, and “” will be displayed throughout the testing.
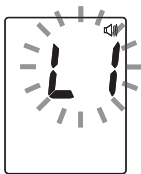
## 5. Choosing a language

Press **M** to select L1/L2/L3.

L1: English.

L2: Spanish.

L3: Chinese.

To confirm your selection, press SET.

## 6. Setting the reminder alarm

Your meter has four reminder alarms.

The meter will display "OFF", "🔔" and "A-1".

If you don't want to set an alarm, press SET to skip this step. Or press

**M** to select "On", then press SET.

With the hour flashing, press **M** to select the correct hour. Press SET.

With the minute flashing, press **M** to select the correct minute. Press SET to confirm and go to the next alarm setting.

**Congratulations! You have completed all settings!**

### NOTE:

- These parameters can **ONLY be changed** in the setting mode.
- If the meter idles for 3 minutes during the setting mode, it will switch off automatically.

# THE FOUR MEASURING MODES



The meter provides you with four modes for measuring, General, AC, PC and QC.

MODES	USE WHEN
General	any time of day without regard to time since last meal
AC	no food intake for at least 8 hours
PC	2 hours after a meal
QC	testing with the control solution


You can switch between the two modes by:

**1.** With the meter turned off, insert a test strip to turn on the meter. The screen will display a flashing blood drop.

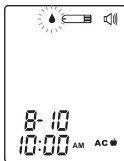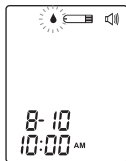**2.** Press **M** to switch between General, AC, PC and QC mode.
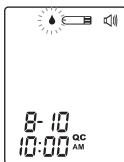
## Control Solution Testing


Our control solutions contain a known amount of glucose that reacts with the test strips, and is used to ensure that the meter and test strips are working together correctly.

Test strips, control solutions, or sterile lancets may not be included in the kit (please check the contents on your product box). They can be purchased separately. Please make sure you have all items prior to testing your blood glucose.

### Test with a control solution when:

- you first receive the meter.
- you do your weekly routine quality testing of the meter and test strips.
- you open a new vial of test strips.
- you suspect the meter or test strips are not working correctly.
- your blood glucose test results are not consistent with how you feel, or if you think the results are not accurate.
- you are practicing the testing procedure.
- you have dropped or think you may have damaged the meter.
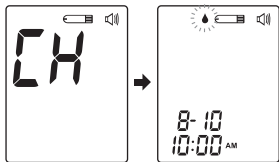
## Performing a Control Solution Test


### 1. Insert the test strip to turn on the meter

Wait for the meter to display “” and “”.

#### Voice Guide

Please relax during measurement. Please apply blood to the strip.

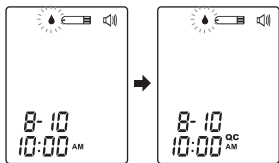
## 2. Press **M** until “QC” is displayed.

With “QC” displayed, the meter will not store your test result in memory. If you press **M** again, the “QC” will disappear and this test is no longer a control solution test.

### **WARNING:**

When doing the control solution test “QC” must be displayed so that the test results will not be stored in the memory. Failure to have the meter in the correct mode will result in having the control solution test results stored with your actual blood glucose test results. This will cause an inaccurate analysis of your stored results.

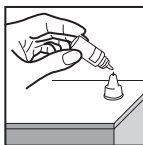
## 3. Apply control solution

Shake the control solution bottle thoroughly before use. Squeeze out a drop onto the tip of the bottle cap and wipe it off using a clean tissue or cloth. Then squeeze out another drop onto the tip of the bottle cap.

Hold the meter so that the strip is pointing downward and move the absorbent hole of the test strip to touch the drop. Once the confirmation window fills completely, the meter will begin counting down.

To avoid contaminating the control solution, do not directly apply control solution onto a strip.
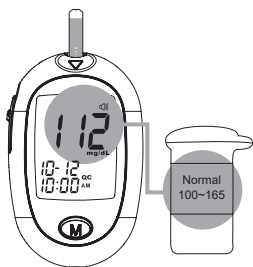
#### 4. Read and compare the result

After counting down to 0, the control solution test result will appear on the display. Compare this result with the range printed on the test strip vial. It should fall within the range. If not, please read the instructions again and repeat the control solution test.

#### Voice Guide

Your blood glucose is (number) milligrams per deciliter.

#### Out-of-range results

If you continue to have test results fall outside the range printed on the test strip vial, the meter and strips may not be working properly. Do **NOT** test your blood using this meter. Go to [www.foracare.com](http://www.foracare.com), or call your customer service for help.

#### NOTE:

- The control solution range printed on the test strip vial is for control solution use only. It is not a recommended range for your blood glucose level.
- See the **Maintenance** section for important information about your control solutions.

# TESTING WITH BLOOD SAMPLE

## Preparing the Lancing Device for Blood Testing

Please follow the instructions in the lancing device insert for collecting a blood sample.

### **WARNING:**

To reduce the chance of infection:

- Never share a lancet or the lancing device.
- Always use a new, sterile lancet. Lancets are for single use only.
- Avoid getting hand lotion, oils, dirt, or debris in or on the lancets or the lancing device.
- Wash and dry your hands thoroughly before and after handling the meter, lancing device and test strips to prevent infection. For more information, please refer to the “Maintenance” section.
- If the meter is being operated by a person who is providing testing assistance to the user, the meter and lancing device should be decontaminated prior to use by that person.

**Sharing the lancing device and lancets may increase the risk of spreading infectious diseases. Lancing device must not be used on more than one person.**

## Preparing the Puncture Site


Rubbing the puncture site before blood extraction has a significant influence on the accuracy of your test results.

Blood from an alternate site that has not been rubbed exhibits a measurably different glucose concentration than blood from the finger. When the puncture site was rubbed prior to blood extraction, the difference was significantly reduced.

**Please follow the suggestions below before obtaining a drop of blood:**


**Wash and dry your hands before starting.**

- Select the puncture site either at fingertips or another body part (please see section “Alternative Site Testing” (AST) on how to select the appropriate sites).
- Rub the puncture site for about 20 seconds before penetration.
- Clean the puncture site using cotton moistened with 70% alcohol and **let it air dry.**
- Use a clear cap while setting up the lancing device.

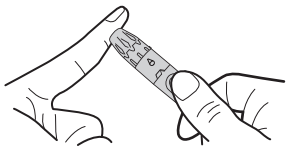- **Fingertip testing**

Press the tip of the lancing device firmly against the lower side of your fingertip. Press the release button to puncture your finger, a click indicates that the puncture is complete.

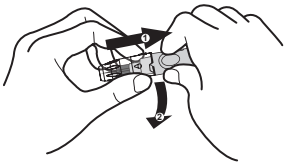
- **Blood from sites other than the fingertip**

Replace the lancing device cap with the clear cap for AST. Pull the cocking control back until it clicks. When lancing the forearm, upper arm, or hand, avoid lancing the areas with obvious veins to avert excessive bleeding.

**NOTE:**

- Choose a different spot each time you test. Repeated punctures at the same spot may cause soreness and calluses.
- Please consult your health care professional before you begin AST.
- It is recommended that you discard the first drop of blood as it might contain tissue fluid, which may affect the test result.

## Alternative Site Testing
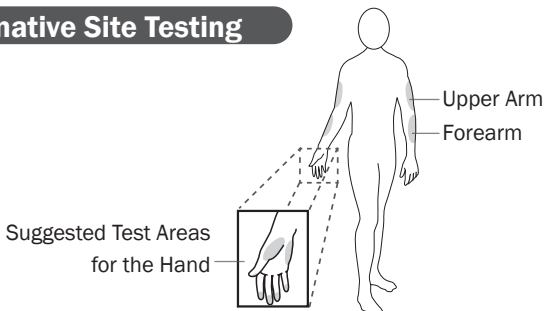
### Important:

**There are limitations with Alternative Site Testing (AST).**

- **AST results should not be used for CGM calibration.**
- **AST results should not be used for insulin dosing calculations.**

**Please consult your health care professional before you perform AST.**

### What is AST?

Alternative site testing (AST) means that people use parts of the body other than the fingertips to check their blood glucose levels. This system provides you to test on the palm, the forearm or the upper arm with the equivalent results to fingertip testing.

### What is the advantage?

Fingertips feel pain more readily because they are full of nerve endings (receptors). At other body sites, since nerve endings are not so condensed, you will not feel as much pain as with your fingertips.

### When to use AST?

Food, medication, illness, stress and exercise can affect blood glucose levels. Capillary blood at the fingertip reflects these changes faster than capillary blood at other sites. Thus, when testing blood glucose during or immediately after a meal, physical exercise, or any other event, **take a blood sample from your finger only.**


We strongly recommend that you perform AST **ONLY** at the following times:

- In a pre-meal or fasting state (more than 2 hours since the last meal).
- Two hours or more after taking insulin.
- Two hours or more after exercise.

Do **NOT** use AST if:

- You think your blood glucose is low.
- You are unaware of hypoglycemia.
- You are testing for hyperglycemia.
- Your AST results do not match the way you feel.
- Your routine glucose results often fluctuate.
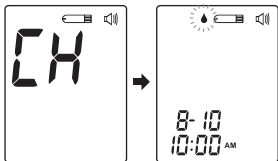
## Performing a Blood Glucose Test

### 1. Insert the test strip to turn on the meter




Wait for the meter to display “” and flashing “”.

#### Voice Guide

Please relax during measurement. Please apply blood to the strip.

### 2. Select the appropriate measuring mode by pressing .

For selecting the measurement mode, please refer to the "THE FOUR MEASURING MODES".
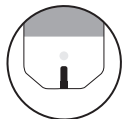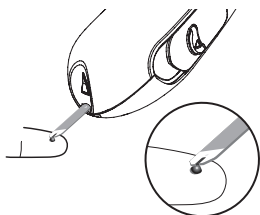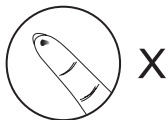
### 3. Obtaining a blood sample

Use the pre-set lancing device to puncture the chosen site. After lancing wipe away the first drop of blood with a clean cotton swab. Gently squeeze the punctured area to obtain another drop of blood. Be careful **NOT** to smear the blood sample.

The volume of blood sample must be at least **0.5** microliter ( $\mu\text{L}$ ) of volume.

( • actual size)

### 4. Apply the sample

Hold the meter with the strip pointing downward and the absorbent hole of the test strip touching the blood drop. Blood will be drawn in to the strip. Once the confirmation window is completely filled, meter begins counting down. Do not remove the strip from the blood drop until the meter beeps.

#### NOTE:

- Do not press the punctured site against the test strip or try to smear the blood.
- If you do not apply a blood sample to the test strip within 3 minutes, the meter will automatically turn off. You must remove and reinsert the test strip to start a new test.
- The confirmation window should be filled with blood before the meter begins to count down. **NEVER** try to add more blood to the test strip after the drop of blood has moved away. **Discard the used test strip and retest with a new one.**
- If you have trouble filling the confirmation window, please contact your health care professional or our customer care service for assistance.


## 5. Read your result

The result of your blood glucose test will appear after the meter counts down to 0. This blood glucose result will automatically be stored in the memory.

### Voice Guide

Your blood glucose is (number) milligrams per deciliter.
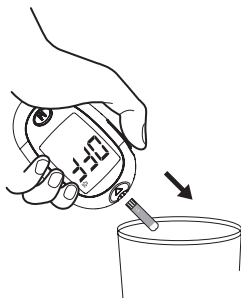
## 6. Eject the used test strip and remove the lancet

To eject the test strip, point the strip into a disposal container for contaminated items. Slide the Test Strip Ejector forward. The meter will automatically turn itself off after the test strip is ejected.

### **Always use with caution when removing the lancet.**

Take the lancet out carefully. Place the disk on a hard surface and push the exposed tip into the protective disk.

**When removing the lancet, follow the instructions from the lancing device insert.**


### **WARNING:**

The used lancet and test strip may be biohazardous. Please follow the proper procedure for discarding biohazardous waste in your area.


# METER MEMORY

The meter stores up to 450 most recent blood glucose test results with date and time stamp in its memory. To enter the meter memory function, **start with the meter turned off.**
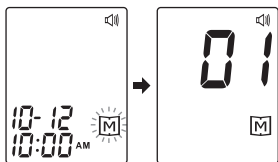
## Reviewing Test Results

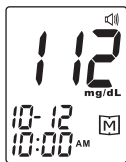
**1. Press and release**  "  " will appear on the display. Press  again, and the first reading you see is the last blood glucose result along with date, time and the measurement mode.

### Voice Guide

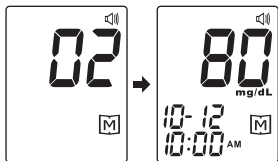Your blood glucose is (number) milligrams per deciliter.

**2. Press**  **again** to recall the next test results stored in the meter each time you press. **After the last test result, press**  **again and the meter will turn off.**


### Voice Guide

Your blood glucose is (number) milligrams per deciliter.

## Reviewing Blood Glucose Day Average Results

**1. Press and release **M**.** When “**M**” appears on the display, keep pressing **M** for 3 seconds until the flashing “**DAY AVG**” appears. Release **M** and then your 7-day average result measured in general mode will appear on the display.

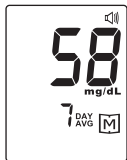**2. Press **M** to review** 14-, 21-, 28-, 60- and 90- day average results stored in memory.

**3. Exit the meter memory.** Press **M** after the 90-day average is displayed and the meter will turn off.

### NOTE:

- Any time you wish to exit the memory, hold press **M** down for 5 seconds or leave the meter idle for 3 minutes. The meter will turn off automatically.
- Control solution results are **NOT** included in the day average.
- If there are no records in memory, “--” displays when you recall the test results or review the average result.


# Downloading Results onto a Computer

## Data Transmission Via RS232 Cable

You can use the meter with FORA 4-pin RS232 cable and the Health Care System Software to view your test results on your personal computer. To learn more about the Health Care System Software or to obtain an RS232 cable, please visit [www.foracare.com](http://www.foracare.com), or call our customer care service number for further assistance.

### 1. Obtaining the required cable and installing the software

To download the Health Care System Software, please visit the ForaCare, Inc. website at <http://www.foracare.com/support.html> and select the “Software & App” to download the “Health Care Software System” software. Follow the online instructions to install the software to your computer.
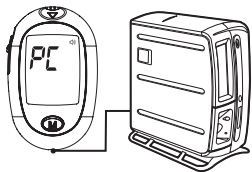
### 2. Connecting to a personal computer

Connect the cable to a RS232 port on your computer. With the meter turned off, connect the other end of the RS232 cable to the meter data port. “PC” will appear on the meter display, indicating that the meter is in communication mode.

### 3. Data transmission

To transmit data, follow the instructions provided with the software. Test results will be transmitted with date and time stamp. After the results have been transmitted, remove the cable, meter will automatically turn off.

#### **WARNING:**

While the meter is connected your computer, it will not be able to perform a blood glucose test.

# MAINTENANCE

## Battery

Your meter comes with two 1.5V AAA size alkaline batteries.

### Low Battery Signal

The meter will display the message below to alert you when the meter power is getting low.


The “” appears with **E-b** and **low**:

The power is not enough to perform testing. Please change the batteries immediately.

### Replacing the Battery

To replace the batteries, make sure that the meter is turned off.

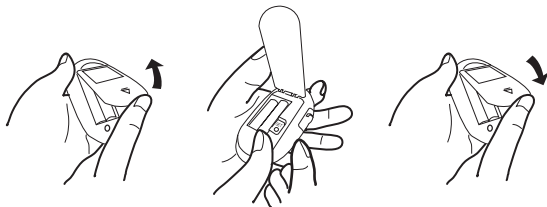1. Remove the battery cover.
2. Remove the old batteries and replace with two 1.5V AAA size alkaline batteries.
3. Close the battery cover, paying close attention to the polarity of the batteries. (If the batteries are inserted correctly, you will hear a “beep”).

#### NOTE:

- Replacing the batteries does not affect the test results stored in the memory.

- As with all small batteries, these batteries should be kept away from children. If swallowed, promptly seek medical assistance.
- Batteries might leak chemicals if unused for a long time. Remove the batteries if you are not going to use the device for an extended period (i.e., 3 months or more).
- Properly dispose of the batteries according to your local environmental regulations.

## Caring for Your Meter

Wash and dry your hands thoroughly before use or handling of the meter and strips.

### What is Cleaning and Disinfection?

Cleaning and disinfection are different. Cleaning is the process of removing dirt (e.g. food debris, grease, dust), disinfection is the process of killing germs (e.g. bacteria and viruses).

### When to clean and disinfect the meter?

Clean the meter when you see any dirt on it. You should disinfect the meter at least once a week to prevent infection.

### How to clean and disinfect the meter?

The meter must be cleaned prior to the disinfection. Use one disinfecting wipe to clean exposed surfaces of the meter thoroughly and remove any visible dirt, blood or any other body fluid with the wipe. Use a second wipe to disinfect the meter. **Do NOT use organic solvents to clean the meter.**

We recommend using the product listed below to clean and disinfect your meter. The recommended product contains isopropyl alcohol as its active ingredient. It has been shown to be safe when used on the FORA


Premium V10 blood glucose meter.

- ▶ Micro-Kill+™ (Micro-Kill Plus™) by Medline  
(EPA Reg. No. 59894-10-37549)

To obtain disinfecting wipes and other information, please contact Medline at 1-800-MEDLINE (1-800-633-5463) or visit [www.medline.com](http://www.medline.com).

## Disinfecting Procedures

1. Take out one disinfecting wipe from the package and squeeze out any excess liquid in order to prevent damage to the meter.

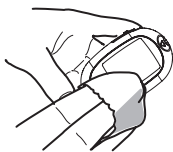2. Wipe all exterior surface of the meter including the display and buttons. Hold the meter with the test strip slot pointing down and wipe the area around the test slot, be careful not to allow excess liquid to get inside. Keep the meter moist with disinfection solution contained in the wipe for a minimum of 2 minutes for Micro-Kill+™ wipes. Follow the instructions on the package label of disinfecting wipe. Use two or more wipes if necessary.

3. Remove the wipe. Allow the meter surface to dry completely.

4. Discard the used wipes and never reuse them.

Wash your hands thoroughly with soap and water after handling the monitor, lancing device and test strips to avoid contamination.

Each cleaning and disinfection cycle includes a pre-cleaning step with one wipe and a disinfection step with a second wipe.

Improper system cleaning and disinfection may result in meter malfunction. If you have a question, please contact customer service at 1-888-307-8188 for assistance.

This device has been validated to withstand 5,000 cleaning and disinfection cycles using the recommended disinfecting wipe/towelette. The tested number of cycles is estimated by 1,000 tests per year on each meter and around 2 cleaning and disinfection cycles per meter per day for 5 years, the expected life of device.

Stop using the meter if you see any signs of deterioration such as the LCD display is cracked or cloudy, the buttons no longer function, or the outer casing is cracked. Please contact customer service at 1-888-307-8188 for a replacement meter if any of the signs of deterioration are noticed.

The meter should be replaced after the validated number of disinfection cycles or the warranty period, whichever comes first.

**NOTE:**

- Do **NOT** clean or disinfect the meter while performing a test.
- Please follow the instructions on the package label of Micro-Kill+™ disinfecting wipes for safe use of the wipes.
- If the meter is being operated by a second person, the meter and lancing device should be decontaminated prior to handling by an assistant.
- Do **NOT** allow cleaning or disinfecting solution to get in the test slot, battery compartment, or strip-ejection button.
- If you do get moisture in the test strip slot, wipe it away with the corner of a tissue.
- Always dry the meter thoroughly before using it.
- Do not spray the meter directly with cleaning solutions especially those containing water (i.e. soapy water), as this could cause the solution to enter the case and damage the electronic components or circuitry.

## Meter Storage

- Storage conditions: -4 °F to 140 °F (-20 °C to 60 °C), below 95% relative humidity.
- Always store or transport the meter in its original storage case.
- Avoid dropping and heavy impact.
- Avoid direct sunlight and high humidity.

## Cleaning and Disinfection of Your Lancing Device

For home users, the lancing device is reusable. Please disinfect your lancing device regularly.

- To clean the lancing device, use a disinfecting wipe to clean the lancing device.
- Disinfect the cap by placing it in 70% alcohol for 10 minutes and then allow it to air dry.
- Do NOT place the device in a dishwasher or use detergents.

## Caring for Your Test Strips

- Storage conditions: 39.2 °F to 104 °F (4 °C to 40 °C), below 85% relative humidity. Do **NOT** freeze.
- Store your test strips in their original vial only. Do not transfer to another container.
- Store test strip packages in a cool dry place. Keep away from direct sunlight and heat.
- After removing a test strip from the vial, immediately close the vial cap tightly.
- Touch the test strip with clean and dry hands.
- Use each test strip immediately after removing it from the vial.
- Write the opening date on the vial label when you first opened it. Discard remaining test strips after 3 months of first opening.

- Do not use test strips beyond the expiration date. This may cause inaccurate results.
- Do not bend, cut, or alter a test strip in any way.
- Keep the strip vial away from children since the cap and the test strip may be a choking hazard. If swallowed, promptly see a doctor for help.

For further information, please refer to the test strip package insert.

## Important Control Solution Information

- Use only our control solutions with your meter.
- Do not use the control solution beyond the expiration date or 3 months after first opening. Write the opening date on the control solution bottle and discard the remaining solution after 3 months of first opening.
- It is recommended that the control solution test be done at room temperature (68°F to 77°F / 20°C to 25°C). Make sure your control solution, meter, and test strips are at this specified temperature range before testing.
- Shake the bottle before use, discard the first drop of control solution, and wipe off the dispenser tip to ensure a pure sample and an accurate result.
- Store the control solution tightly closed at temperatures between 35.6°F and 86°F (2°C and 30°C). Do **NOT** freeze.

# SYSTEM TROUBLESHOOTING

If you follow the recommended action but the problem persists, or error messages other than the ones below appear, please call your local customer service. Do not attempt to repair the meter yourself and never try to disassemble the meter under any circumstances.

## Result Readings

MESSAGE	WHAT IT MEANS
 The image shows a digital display with 'Lo' in large characters, a small 'dL' in the top right, and the time '10:12' and '10:00 AM' at the bottom.	Appears when your blood glucose test result is below the lower measurement limit, which is less than 20 mg/dL.
 The image shows a digital display with '388' in large characters, 'mg/dL' and 'KETONES' in smaller characters below it, and the time '10:12' and '10:00 AM' at the bottom.	Appears when your result is equal to or higher than 240 mg/dL. This indicates the possibility of ketone accumulation for type 1 diabetes. Please seek medical assistance immediately.
 The image shows a digital display with 'Hi' in large characters, a small 'dL' in the top right, and the time '10:12' and '10:00 AM' at the bottom.	Appears when your blood glucose test result is higher than the upper limit of measurement, which is 600 mg/dL.

## Error Messages

MESSAGE	WHAT DOES METER SAY?	WHAT IT MEANS	WHAT TO DO
	Only Error Beep	Problem with the meter.	Review the instructions and repeat the test. If the meter still does not work, contact the local distributor for assistance.
	Only Error Beep		
	Only Error Beep		
	Only Error Beep		
	Only Error Beep	Problem with the meter.	Repeat the test with a new test strip. If the meter still does not work, please contact the customer service for assistance.
	Only Error Beep	Appears when the batteries cannot provide enough power for a test.	Replace the batteries immediately.

MESSAGE	WHAT DOES METER SAY?	WHAT IT MEANS	WHAT TO DO
	Only Error Beep	Appears when test strip is removed while counting down.	Repeat the test with a new test strip.
	Only Error Beep	Appears when a used test strip is inserted.	Repeat with a new test strip.
	Only Error Beep	Appears when ambient temperature is below system operation range.	System operation range is 50 °F to 104 °F (10 °C to 40 °C). Repeat the test after the meter and test strip are in the above temperature range.
	Only Error Beep	Appears when ambient temperature is above system operation range.	

## Troubleshooting

1. If the meter does not display a message after inserting a test strip:

POSSIBLE CAUSE	WHAT TO DO
Batteries exhausted.	Replace the batteries.
Test strip inserted upside down or incompletely.	Insert the test strip with contact bars end first and facing up.
Defective meter or test strips.	Please contact customer service.

2. If the test does not start counting down after applying the sample:

POSSIBLE CAUSE	WHAT TO DO
Defective test strip.	Repeat the test with a new test strip.
Sample applied after automatic switch-off (3 minutes after last user action).	Repeat the test with a new test strip. Apply sample only when flashing " ♠ " appears on the display.
Defective meter.	Please contact customer service.

3. If the control solution testing result is out of range:

POSSIBLE CAUSE	WHAT TO DO
Error in performing the test.	Read instructions thoroughly and repeat the test again.
Control solution vial was poorly shaken.	Shake the control solution vigorously and repeat the test again.
Expired or contaminated control solution.	Check the expiration date of the control solution.
Control solution that is too warm or too cold.	Control solution, meter, and test strips should be at room temperature (20 °C to 25 °C/68 °F to 77 °F) before testing.
Defective test strip.	Repeat the test with a new test strip.
Meter malfunction.	Please contact customer service.

# DETAILED INFORMATION

The meter provides you with plasma equivalent results.

Time of day	Normal plasma glucose range for people <b>without</b> diabetes (mg/dL)
Fasting and before meal	< 100 mg/dL (5.6 mmol/L)
2 hours after meals	< 140 mg/dL (7.8 mmol/L)

Source: American Diabetes Association. Standards of medical care in diabetes-2016; 39 (supp.1 Diabetes Care): S16.

**Please consult your doctor to determine a target range that works best for you.**

# SPECIFICATIONS

**Model No.:** Premium V10

**Dimension & Weight:** 88.67 x 55.61 x 23.07 mm, 48.8 g

**Power Source:** Two 1.5V AAA alkaline batteries

**Display:** LCD

**Memory:** 450 measurement results with respective date and time

**External output:** RS232

Auto electrode insertion detection

Auto reaction time count-down

Auto switch-off after 3 minutes without action

Temperature Warning

**Operating Condition:** 50 °F to 104 °F (10 °C to 40 °C), below 85% R.H.  
(non-condensing)

**Meter Storage/Transportation Conditions:**

-4 °F to 140 °F (-20 °C to 60 °C), below 95% R.H.

**Strip Storage/Transportation Conditions:**

39.2 °F to 104 °F (4 °C to 40 °C), below 85% R.H.

**Measurement Units:** mg/dL

**Measurement Range:** 20 to 600 mg/dL

This device has been tested to meet the electrical and safety requirements of: IEC/EN 61010-1, IEC/EN 61010-2-101, EN 61326-1, IEC/EN 61326-2-6.

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**Distributed by ForaCare, Inc. / Distribuido por ForaCare, Inc.**